

DANIEL POCCI
ITALMATCH CHEMICALS
1100 STATE ROAD
ASHTABULA OH 44004
UNITED STATES

Phone: 216-767-4001
Email: D.POCCI@italmatch.com

Work Order: 23-07762
Revision: 1
Report Date: 11-Sep-23
Order Date: 21-Aug-23
P.O.: 70318334

Release:
Approved By: Lanny Huynh
Title: Manager, Organic Analysis

BALAZS™ TEST RESULTS

If you have any questions regarding the results, please call Tiffany Wilkus at (510) 624-4045

The information contained in this electronic or original document is confidential information for the intended recipient. If you are not the intended recipient, or the employee or agent responsible for delivering this document to the intended recipient, you are hereby notified that any disclosure, dissemination, distribution or copying of this communication is strictly prohibited. If you have received this document in error, please immediately notify the sender and destroy this document.

Background and Summary:

On 8/21/2023, one HCl solution (Hydrochloric acid, LOT 081523114) was received from Italmatch Chemicals. The sample was extracted with hexane and the extracted solution was then analyzed by GC-MSD (Gas Chromatography with Mass Selective Detector) for qualitative and semi-quantitative analysis of extractable organic contaminants, and GC-FID (Gas Chromatography with Flame Ionization Detector) for quantitative determination of chlorobenzene.

Sample ID	Sample Mass (g)
LOT 081523114	33.62

Sample Preparation:

Approximately 30mL of sample was extracted with 3mL of hexanes. The extracted solution was concentrated to 1mL and spiked with a known volume of n-Hexadecane internal standard. The extract was analyzed by GC-MS and GC-FID.

A calibration curve for chlorobenzene ranging in concentration from 0.1 - 100 ppb was also run by GC-FID.

GC-FID Instrumentation:

The system used for the analysis was an Agilent 7890B GC with a FID. The GC-FID was equipped with a non-polar phase (dimethyl polysiloxane) capillary column. The following temperature program was used for the GC: the oven was held at an initial temperature of 35°C for 5 minutes, ramped at a rate of 15°C/minute

These results were obtained by following standard laboratory procedures and are only representative of the sample as received by the laboratory. The liability of AIR LIQUIDE - BALAZS NanoAnalysis ("Balazs") should not exceed the amount paid for this report. In no event shall Balazs be liable for special or consequential damages. Client agrees not to use Balazs's name in reporting results obtained from tests performed by Balazs written consent as to such use. Report shall not be reproduced except in full, without the written permission of Balazs.

to 280°C and held at the final temperature for 20 minutes. The injection port was operated in splitless mode at a temperature of 250°C. The sample injection volume was 1.0 µL.

GC-FID Results:

Quantitative results: Quantitative amounts of the chlorobenzene were calculated using the calibration curve. Chlorobenzene was not detected at or above 10 ppbv. No other impurities were detected by this method.

GC-MS Instrumentation:

The system used for the analysis was an Agilent 7890B GC with a 5977 quadrupole MSD. The GC was equipped with a non-polar dimethyl polysiloxane phase capillary column. The following temperature program was used for the GC: the oven was held at an initial temperature of 35°C for 5 minutes, ramped at a rate of 15°C/minute to 280°C and held at the final temperature for 25 minutes. The injection port was operated in splitless and a temperature of 230°C. The sample injection volume was 1.0 µL.

Each compound passed down the GC column at a characteristic rate. As each compound exited the gas chromatograph, it entered the MSD where it was ionized using electron impact ionization (70 eV). The MSD collected a full mass spectrum (10-700 amu) approximately once per second.

GC-MS Results:

Identification of each compound detected was first attempted by searching a Wiley/NIST library of nearly 900,000 mass spectra. In cases where no matches were found, the analyst interpreted the mass spectra to provide the best estimate of the most probable compound or class of compounds.

Labeled chromatograms are included in the report.



Air Liquide Electronics U.S. LP
46409 Landing Parkway, Fremont CA 94538
Telephone (510) 624-4000 Fax (510) 657-2292

DANIEL POCCHI
ITALMATCH CHEMICALS
P.O.: 70318334
Release:

Report Date: 11-Sep-23
Order Date: 21-Aug-23
Work Order: 23-07762
Approved By: Lanny Huynh
Title: Manager, Organic

Sample: 1417890 [LOT 081523114]

Site ID: HCl

Component	Detection Limit	Result Value	Units
G0111-GCL-N-R-C Solvent (Mixture) Normalized Area % Assay by GC- FID			
Chlorobenzene	10	*	ppb

* = Analysis revealed that the analyte was not found at or about the detection limit

These results were obtained by following standard laboratory procedures and are only representative of the sample as received by the laboratory. The liability of AIR LIQUIDE - BALAZS NanoAnalysis ("Balazs") should not exceed the amount paid for this report. In no event shall Balazs be liable for special of consequential damages. Client agrees not to use Balaz's name in reporting results obtained from tests performed by Balazs written consent as to such use. Report shall not be reproduced except in full, without the written permission of Balazs.

Figure 1: Chromatogram (Overlays-Full Scale)


Sample ID: Hexane Extraction Blank, HCl Lot 081523114

Method of Analysis: Liquid Injection-GC-MS (Gas Chromatography-Mass Spectrometry)

Sampling Conditions: Direct injection of hexane extract

Internal Standard: n-Hexadecane

Balazs Work Order: 23-07762

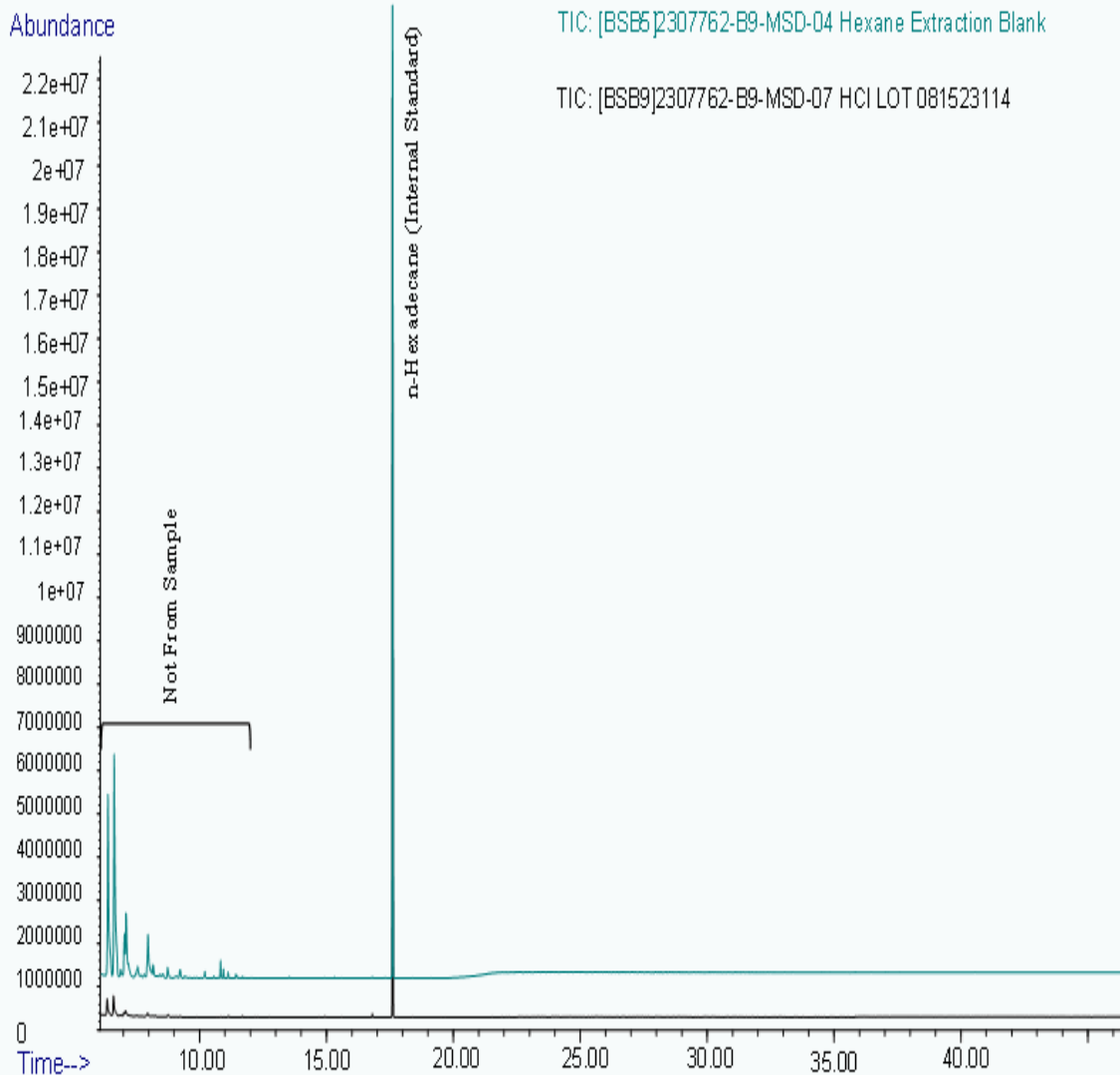

These results were obtained by following standard laboratory procedures and are only representative of the sample as received by the laboratory. The liability of AIR LIQUIDE - BALAZS NanoAnalysis ("Balazs") should not exceed the amount paid for this report. In no event shall Balazs be liable for special or consequential damages. Client agrees not to use Balazs's name in reporting results obtained from tests performed by Balazs written consent as to such use. Report shall not be reproduced except in full, without the written permission of Balazs.

Figure 2: Chromatogram (Overlays - Expanded/zoomed into impurities)


Sample ID: Hexane Extraction Blank, HCl Lot 081523114

Method of Analysis: Liquid Injection-GC-MS (Gas Chromatography-Mass Spectrometry)

Sampling Conditions: Direct injection of hexane extract

Internal Standard: n-Hexadecane

Balazs Work Order: 23-07762

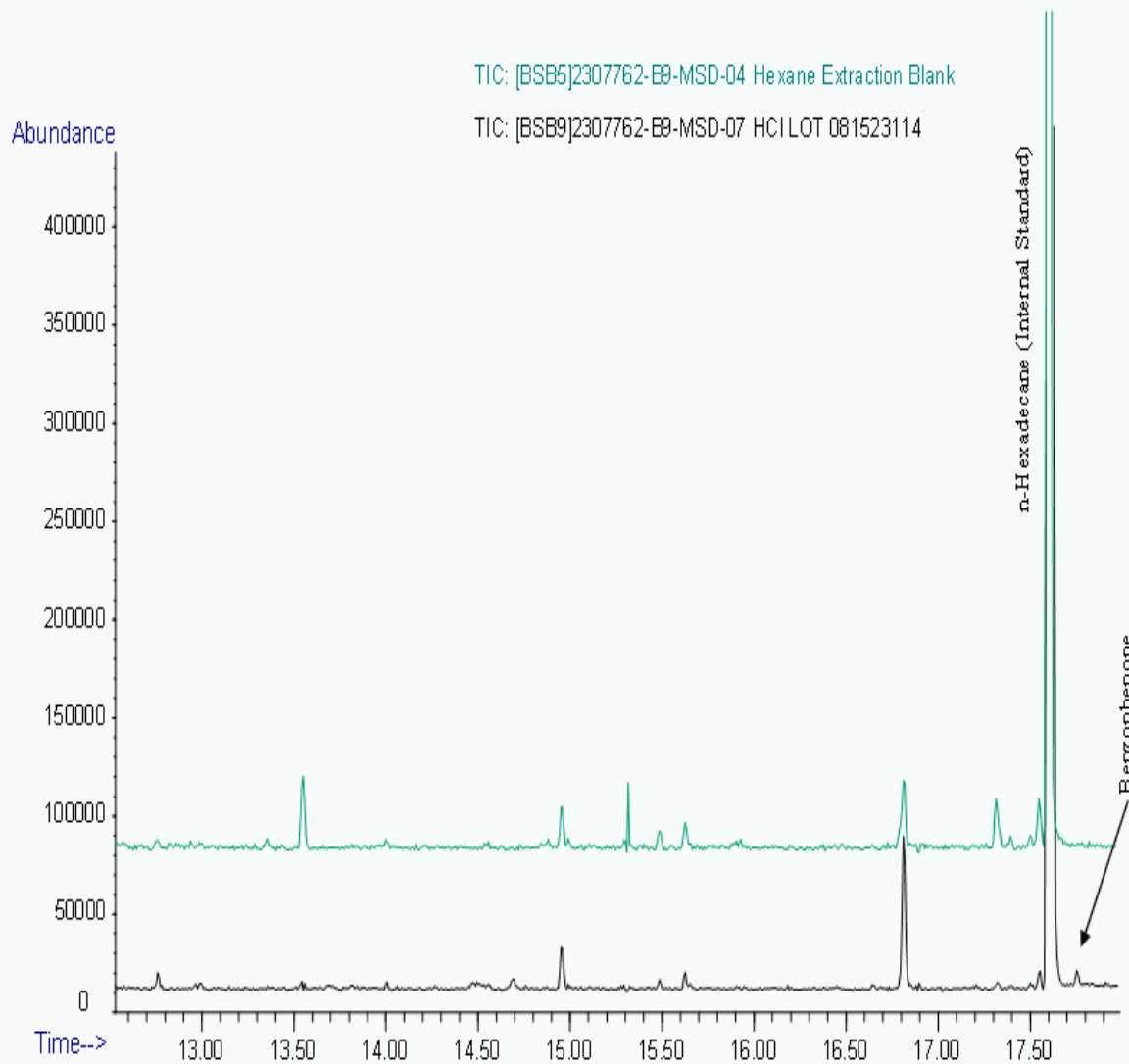

These results were obtained by following standard laboratory procedures and are only representative of the sample as received by the laboratory. The liability of AIR LIQUIDE - BALAZS NanoAnalysis ("Balazs") should not exceed the amount paid for this report. In no event shall Balazs be liable for special or consequential damages. Client agrees not to use Balazs's name in reporting results obtained from tests performed by Balazs written consent as to such use. Report shall not be reproduced except in full, without the written permission of Balazs.

Figure 3: Chromatogram (Overlays - Expanded into impurities)


Sample ID: Hexane Extraction Blank, HCl Lot 081523114

Method of Analysis: Liquid Injection-GC-MS (Gas Chromatography-Mass Spectrometry)

Sampling Conditions: Direct injection of hexane extract

Internal Standard: n-Hexadecane

Balazs Work Order: 23-07762

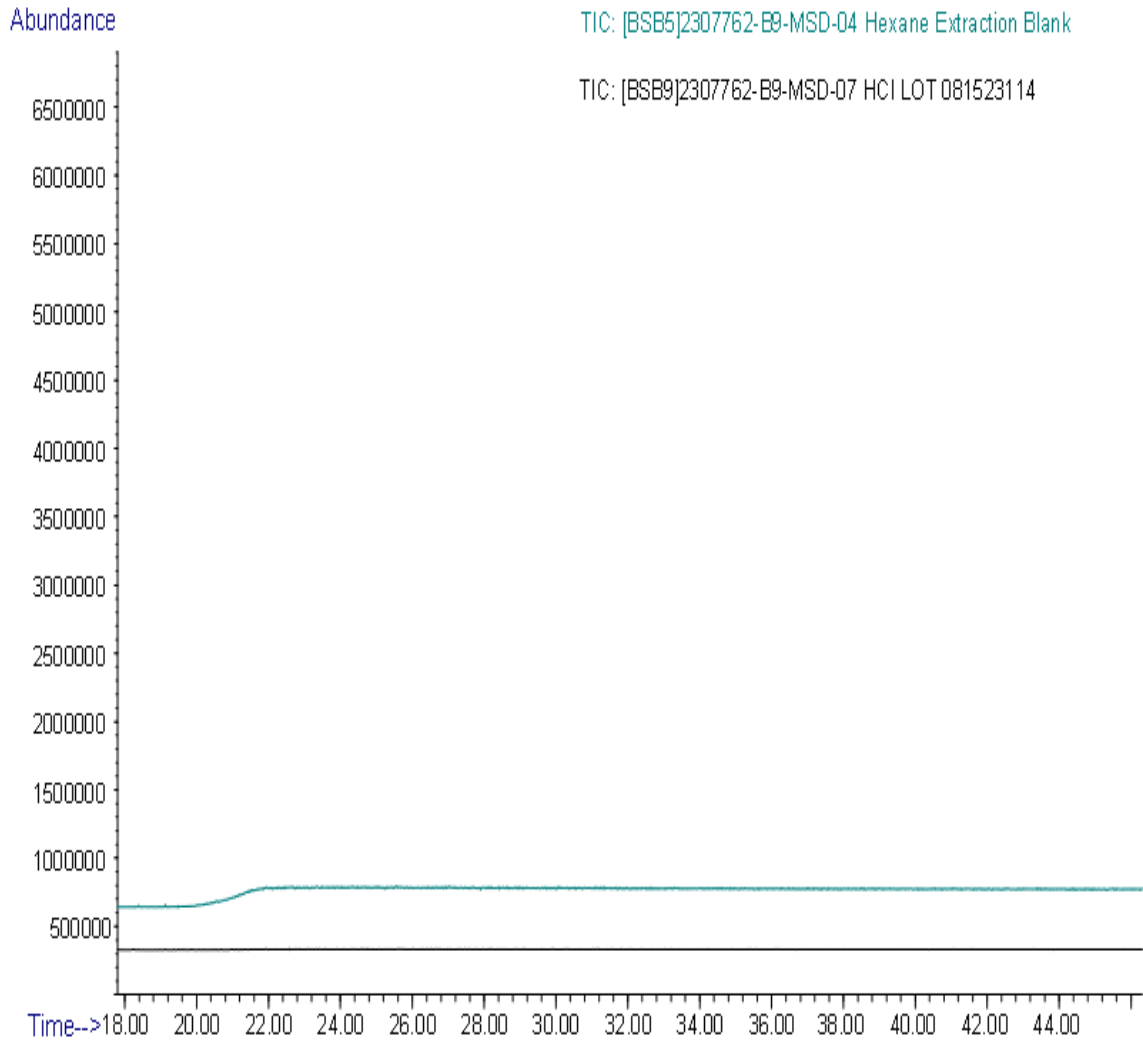

These results were obtained by following standard laboratory procedures and are only representative of the sample as received by the laboratory. The liability of AIR LIQUIDE - BALAZS NanoAnalysis ("Balazs") should not exceed the amount paid for this report. In no event shall Balazs be liable for special or consequential damages. Client agrees not to use Balazs's name in reporting results obtained from tests performed by Balazs written consent as to such use. Report shall not be reproduced except in full, without the written permission of Balazs.



Air Liquide Electronics U.S. LP

46409 Landing Parkway, Fremont CA 94538

Telephone (510) 624-4000 Fax (510) 657-2292

DANIEL POCCHI
ITALMATCH CHEMICALS

P.O.: 70318334

Release:

Report Date: 11-Sep-23

Order Date: 21-Aug-23

Work Order: 23-07762

Approved By: Lanny Huynh

Title: Manager, Organic

These results were obtained by following standard laboratory procedures and are only representative of the sample as received by the laboratory. The liability of AIR LIQUIDE - BALAZS NanoAnalysis ("Balazs") should not exceed the amount paid for this report. In no event shall Balazs be liable for special or consequential damages. Client agrees not to use Balaz's name in reporting results obtained from tests performed by Balazs written consent as to such use. Report shall not be reproduced except in full, without the written permission of Balazs.